

lifelong SAFEWAY

SAFEWAY SAFETY SYRINGE

Latex Free : Pyrogen Free : Non Toxic


Safeway Safety Syringe protect health professionals from “needlestick” hazard because after use the contaminated needle locked in safety cover.

PRODUCT FEATURE & QUALITY:

- Manufactured as per IS : 10654 / ISO 7864.
- Safety Syringe confirms to ISO 7864-1, ISO 594-1, ISO 594-2, ISO 3696.
- Made of plastic polymer inert, in other words, does not react with the medication.
- Sterilized by Ethylene Oxide (ETO).
- Body Transparent - allows perfect visualization of the level of medicine, besides facilitating the detection of possible impurities.
- Smooth plunger - lubrication with high-quality silicone reduces the required pressure during the injection, making it much more comfortable.
- Slider with lock - prevents the separation between the plunger and the syringe body, preventing leaks.
- Nozzle luer-lock.
- **Provided with a safety device that covers the needle after activation, minimizing risk of injury.**




ACTIVATION OF THE SAFETY CLIP:

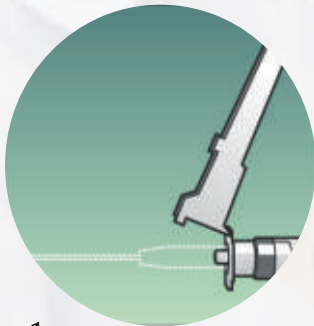


Step 1.

The contaminated needle is connected to the safety Syringe.

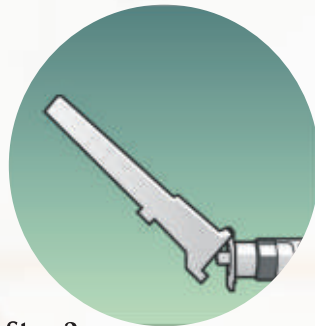

Step 2.

To activate the safety device, push the device on the needle.

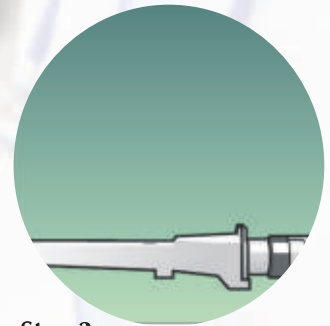

Step 3.

With the activated device, the needle is covered.

Product Range

Available in size 3-5-10-20ml luerlock

Lifelong Meditech Ltd.

Corporate Office : D-1, Soami Nagar (South), New Delhi-110017, India • Tel. : +91-11-46017700-98

Works : Plot No. 23, Sector-5, IMT Manesar, Gurgaon (Haryana)-122050, India • Tel. : +91-124-4335800 • Fax : +91-124-4335887

Website : www.lifelongmed.com • Email : info@lifelongindia.com • Customer Care : lmlcustomercare@lifelongindia.com